


Mobile digital baby scale

seca 334

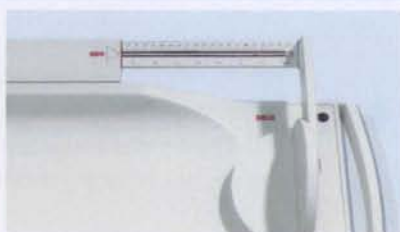

Handy and light for mobile medical use . Comfortable cradle . Large, easy to read LCD display . Reliable results due to automatic damping system . Optional measuring rod seca 232

Because babies move around – and doctors have to be able to move around. The mobile baby scales with special damping system.

Modern technology can only be of use if it passes the tests of everyday life. The digital baby scale seca 334 has all that is necessary to make daily routine in the surgery and hospital so much easier in future: as this scale is not only easy to use but also mobile and durable, it meets all the prerequisites of an ideal scale for paediatrics.



Handy and lightweight for mobile medical use

Weighing only 3.5 kg, this scale is a real lightweight in its class. The integrated side handle makes it easy to transport. Another practical feature is the flat and thus space-saving design, making it easy, for example, to store the scale away in the drawer of a baby's changing table. It can be packed away safely for transport in the optional transport case seca 428.




Easy to operate with generously sized cradle and easy to read LCD display

The baby feels at ease and safe in the extra large, comfortable weighing tray of the seca 334. Larger babies or small children can be weighed in a sitting position. The weight is shown on the brightly illuminated LCD display and the automatic switch-off ensures power-saving operation.



Reliable results due to the automatic damping system

Absolute precision in all situations: fitted with highly sensitive weighing technology and the special seca damping system, exact results are a matter of course for the seca 334. Even when babies are very restless.

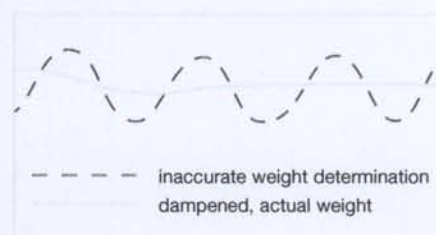

Optional measuring rod seca 232

When fitted with the optional measuring rod seca 232, the baby scale seca 334 becomes a time-saving, combined measuring and weighing station.



Technical data

seca 334

Capacity: 20 kg / 44 lbs

Graduation: 5 g < 10 kg > 10 g /

0.2 oz < 22 lbs > 0.5 oz

Dimensions (WxHxD): 640 x 105 x 300 mm

Weight: 3.5 kg

Functions: Automatic switch-off, TARE, auto-Hold, kg/lbs switch-over

Optional: Measuring rod seca 232, head positioner seca 419, transport case seca 428, AC adapter